


General Laparoscopy Tray and Long Laparoscopy Tray




**Custom Protection,
Ultimate Security**

Custom Protection, Ultimate Security

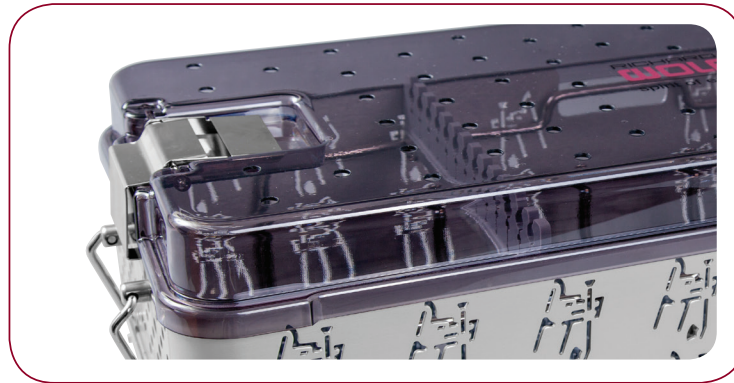


Domed lid accommodates Richard Wolf Eragon style laparoscopy instruments

Cut-outs on the sides of the tray allow for easy transfer of the instrument carrier from the tray

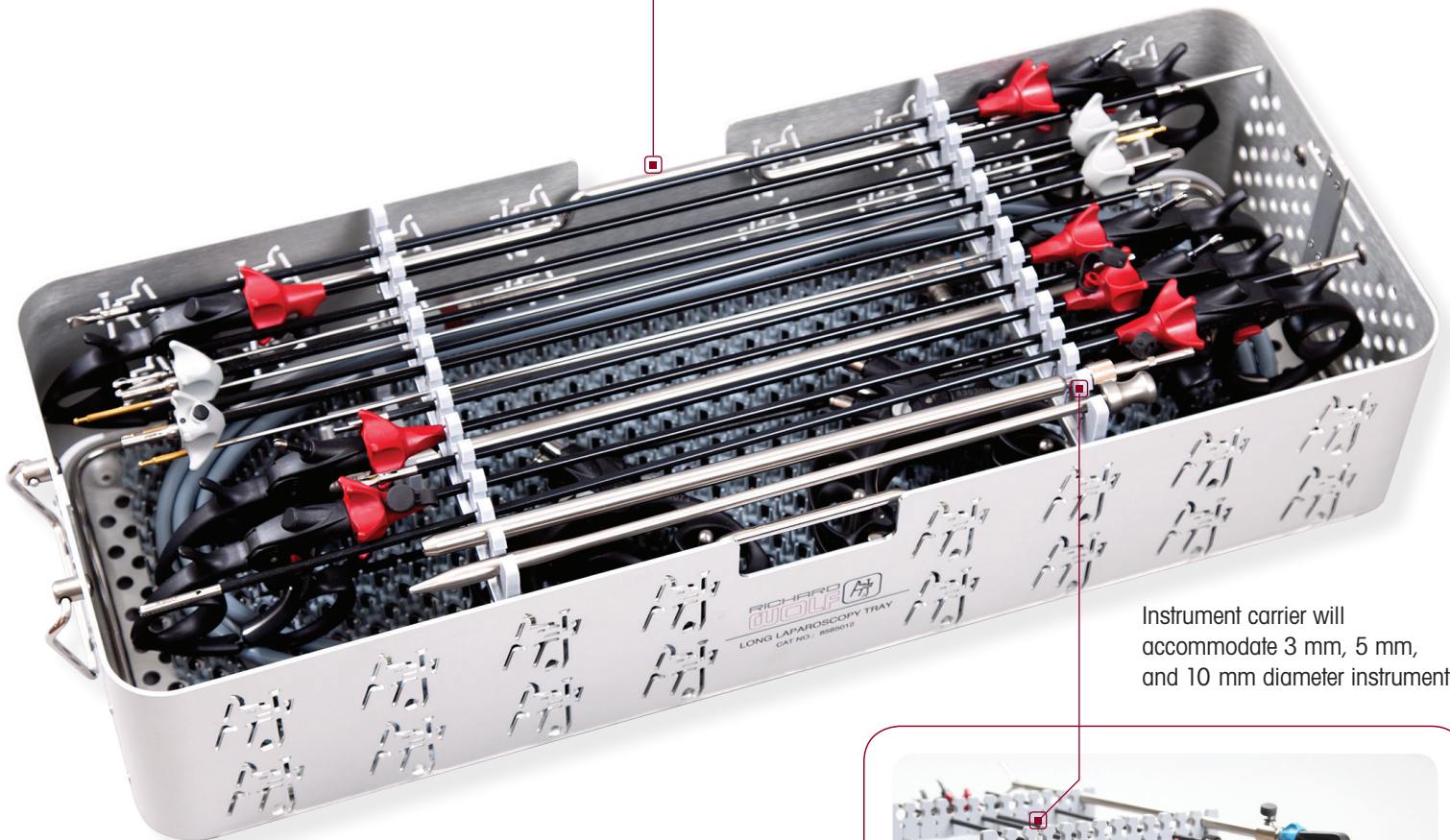


Instrument carrier will accommodate 3 mm, 5 mm, and 10 mm diameter instruments

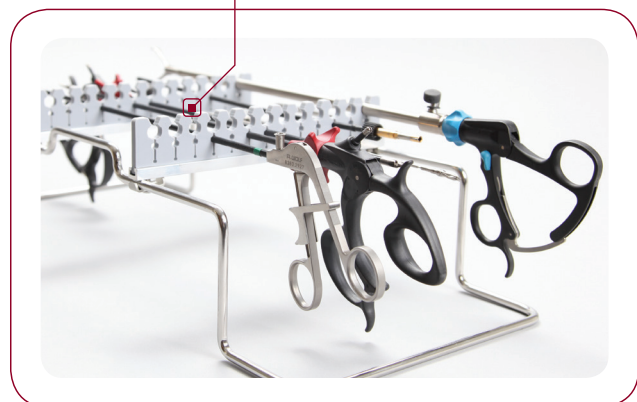

Custom Protection, Ultimate Security

- Aluminum posts hold up to 6 handles used with modular laparoscopy instruments.
- Cut-outs on the sides of the tray allow for easy hand access to remove the instrument carrier from the tray.
- Light cable fits in the perimeter of the tray.
- Deep pin matting secures miscellaneous instruments.
- Brackets to secure laparoscopes

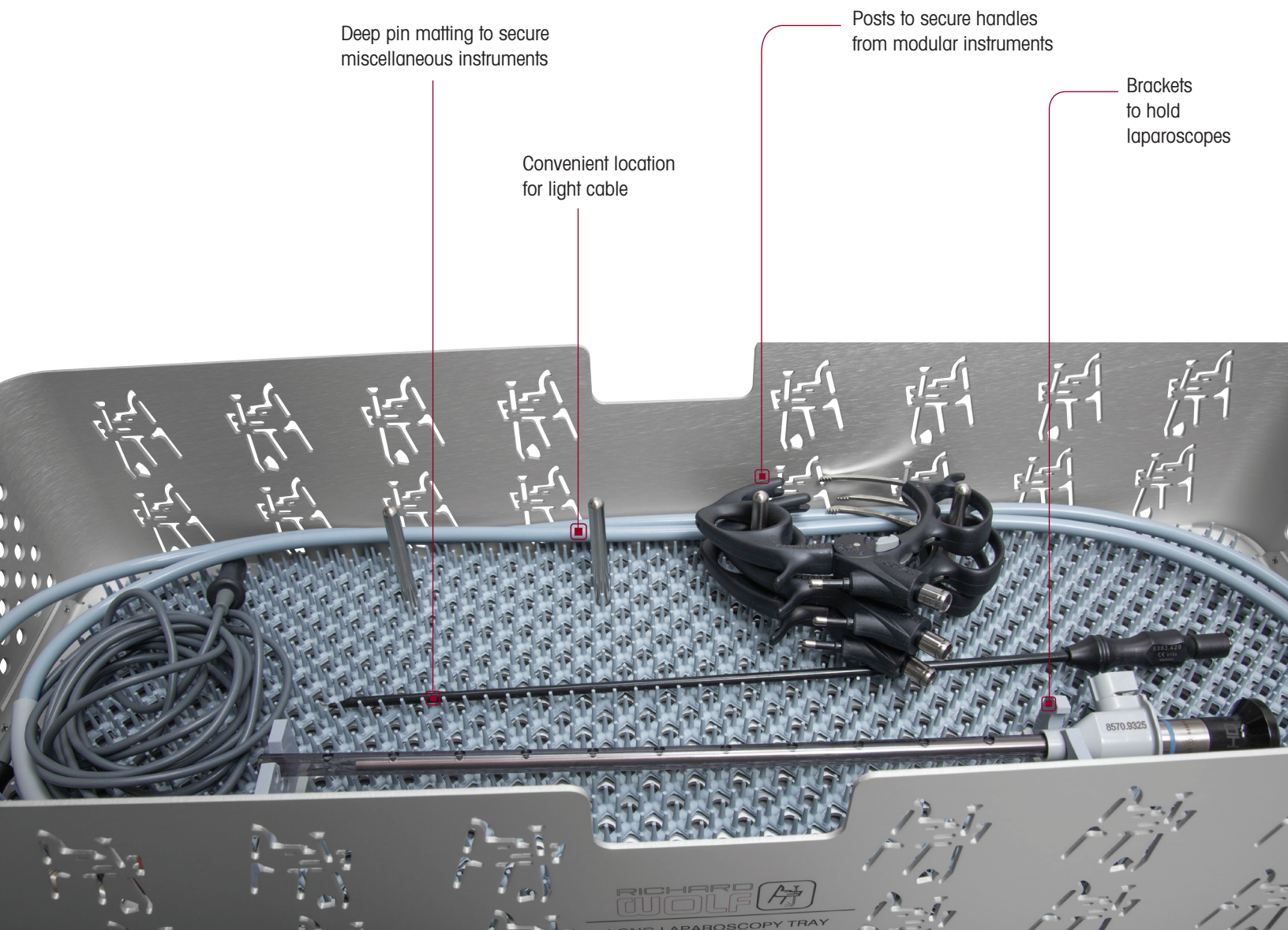

Custom Protection, Ultimate Security

Key Features:

General Laparoscopy Tray



The General Laparoscopy Trays is indicated for packaging and sterilizing the following:

General Laparoscopy Tray

- Up to 13 single-piece or modular standard length (320 mm) laparoscopy instruments
- Fits 3mm, 5mm or 10mm diameter instruments
- 1 or 2 standard length (320 mm) laparoscopes, and a light cable
- Up to 6 modular handles for modular instruments
- Options for miscellaneous cautery tips, electrical cord (monopolar or bipolar)

Long Laparoscopy Tray



The Long Laparoscopy Tray is indicated for packaging and sterilizing the following:

Long Laparoscopy Tray

- Up to 13 single-piece or modular laparoscopy instruments fits either standard length (320 mm) or long length (420 mm – 440 mm)
- Fits 3 mm, 5 mm or 10 mm diameter instruments
- 1 or 2 laparoscopes, either standard length (320 mm), long length (420 mm – 440 mm) or combination of both, and a light cable
- Up to 6 modular handles for modular instruments
- Options for miscellaneous cautery tips, electrical cord (monopolar or bipolar)

Ordering Information

Tray Name	Catalog Number	Dimensions	Sterilization Validation
General Laparoscopy	8585010	20.4" L x 9.6" W x 4.89" D	Prevac Steam
Long Laparoscopy	8585012	24.95" L x 9.6" W x 4.89" D	Prevac Steam