

The Future of Diagnostic Thoracoscopy



Mini Thoracoscopy Set


for Single Port Technique

Mini Thoracoscopy Set

for Single Port Technique

This all-in-one Instrument Set for Single Puncture Thoracoscopy allows minimally invasive diagnostic treatment with one small incision. With the small trocar the procedure is possible under local anesthetics and without painful stress for the patient after surgery.

- Very small Thoracoscopy Set with only 5.5 mm diameter
- Optimum image quality with 50.000 pixels
- Large working channel for 3.5 mm instruments
- Very atraumatic method for the patient - reduced pain
- Suction of fluids possible with instrument inserted
- No conventional rod lenses means lens breakage minimized
- The biopsy forceps provide high quality tissue samples
- The working instruments are very easily changed and are always in the field of view
- Ergonomic design

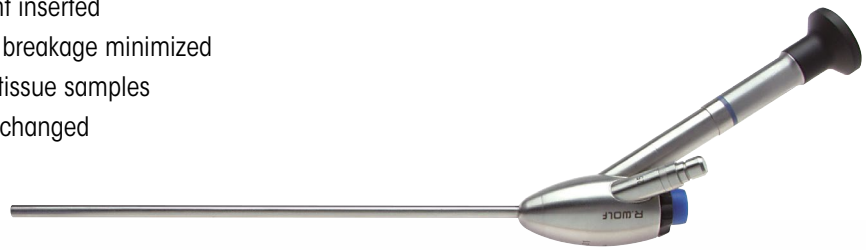

OP-Endoscope bndl "Leipzig Egg"

Comprising:

5.5 mm Operating Laparoscope with 3.5 mm working channel, direction of view 0°, WL 215 mm, silicate image guide (8920.401), sealing membrane, Ø 17 mm (89.103), seal (8920.311), seal cap (15176.020)..... 8920.4011

Flexible Trocar Sleeve, Ø 5.5 mm,

Plastic sheath with thread, distal straight, WL 60 mm 8905.331

also:

Trocar, spherical 8905.3315

Flexible Trocar Sleeve, Ø 5.5 mm,

Plastic, WL 75 mm 8906.051

also:

Trocar,

Tip tapered blunt 8906.151

Double spoon forceps monopolar

Ø 3.5 mm bndl, modular, spoon-shaped, WL 310 mm, insulated
Comprising:

Handle pistol shaped mono (83930074), sheath tube insulated, Ø 3.5 mm, SL 330 mm (8391933), double spoon forceps insert, Ø 3.5 mm (8391216)..... 83912167

Hook Electrode,

monopolar, Ø 3.5 mm, WL 310 mm 8379.452

Coagulation Button Electrode,

monopolar, Ø 3.5 mm, WL 310 mm 8379.462

Probe Rod,

Ø 3.5 mm, graduated, WL 350 mm 8379.672

Suction Tube,

Ø 3 mm, graduated, WL 450 mm 8380.68

Fiber light cable bndl,

Ø 2.5 mm, TL 2.3 m

Comprising:

Fiber light cable (80662523), adapter projector side (8095.07), adapter endoscope side (809509)..... 806625231

Perforated basket for telescope

Holds: Telescope, for machine reprocessing, sterilization, storage and transport 38008

Subject to modifications. Please note that some products may not be available in all countries. Please always take note of the information that is provided with the product. For more detailed information, please consult with your contact partner or directly with Richard Wolf GmbH.

B779en230224

