


## The next Generation of Rigid Bronchoscopes

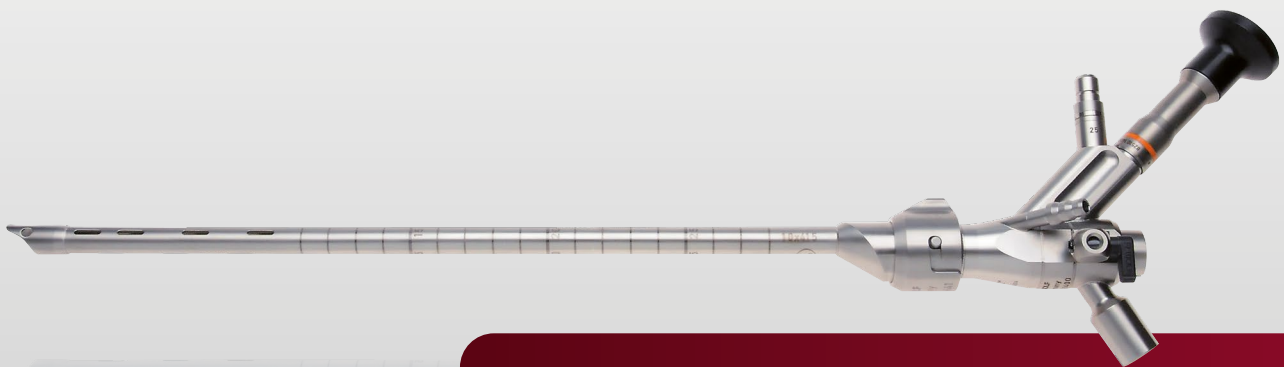

### **TEXAS**

The First Optical  
Fully Integrated Rigid Bronchoscope

## Optical Fully Integrated Rigid Bronchoscope

The **TEXAS Optical Fully Integrated Rigid Bronchoscope** was designed in close cooperation with Dr. Garrett Walsh\* for both diagnostic and therapeutic bronchoscopy applications, such as biopsy, tumor resection, foreign body removal, stent placement and removal. The integrated telescope adds another layer of protection, with separate channels for optic and instruments.



### TEXAS Semiflexible Endoscope

Direction of view 17°, with integrated lens irrigation and protective sheath

(15208.289) .....82510.1400

### TEXAS Bronchoscope Tubes



Size / Inner ø	Outer ø	WL	Type
6 mm	10 mm	305 mm	82520.0641
8 mm	12 mm	310 mm	82520.0841
8 mm	12 mm	310 mm	82520.0842**
10 mm	14 mm	318 mm	82520.1041
12 mm	16 mm	325 mm	82520.1241
14 mm	18 mm	330 mm	82520.1441

- The endoscope is integrated directly in the bronchoscope tube and always offers an optimum view.
- Intubation is made significantly easier as the view is always available.
- The integrated irrigation channel means that the distal lens of the endoscope can be irrigated during an operation if it becomes cloudy due to secretions or blood.
- Easy changeover between different bronchoscope tubes.
- Compatible with all forceps intended for this purpose and for combined application (rigid and flexible working) with the flexible bronchoscope.
- The supplementary range of instruments used, e.g. foreign-body forceps, can be quickly changed and is always in the optical field of view.

### Foreign Body Forceps,

with alligator jaws .....8280.41

**Universal Forceps** .....8280.42

**Biopsy Forceps** .....8280.43

### Grasping Forceps,

for soft foreign bodies .....8280.47

**Laser Guide Tube** ..... 82510.0011

### Nozzle for Jet Ventilation bundle

including luer connector

(15401.071) ..... 8238.5021

**Universal Lock Cap** .....8020.15

### Fiber light cable bndl,

including Fiber light cable ø 2.5 mm TL 2.3 m (80662523), Adapter endoscope side (809509), Adapter projector side (809507)..... 806625231

\* Dr. Garrett Walsh, Professor of Surgery, Dept. of Thoracic and Cardiovascular Surgery, University of Texas MD Anderson Cancer Center, Houston, TX

\*\* without distal window