

Exceptional Image Quality
in a Smaller Design




3.8 Compact Hysteroscope

Slim 3.8 mm outer diameter
for less dilation

3.8 mm Compact Hysteroscope

features increased patient comfort, enhanced durability, and effortless design while maintaining excellent image quality.

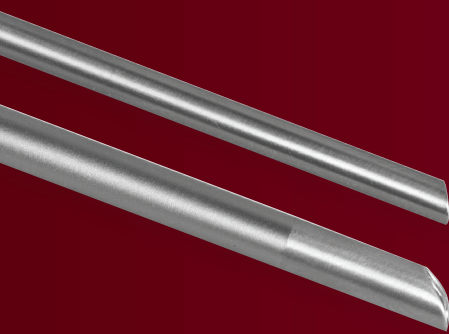

- Approximately 30% thinner than comparable hysteroscopes
- Slim 3.8 mm outer diameter requires less dilation than conventional hysteroscopes
- Rod lens optics system offers superior image quality over flexible and rigid fiber optic scopes
- 5 Fr. working channel accommodates a larger range of instruments for expanded treatment options
- Dedicated inflow and outflow channels maintain distention and visibility throughout the entire procedure

Hysteroscopes				
Catalog Number	Outer Diameter	Viewing Angle	Working Channel	Working Length
8986402	3.8 mm	30°	5 Fr.	217 mm
8986.401	5 mm	20°	5 Fr.	217 mm
HySafe Qui-Connect Modular Forceps, Rotatable				
Catalog Number	Description			
8642.6002	HySafe spoon biopsy forcep, 5 Fr., 340 mm WL			
8642.6312	HySafe punch biopsy forcep, 5 Fr., 340 mm WL			
8642.6502	HySafe Alligator grasper, 5 Fr., 340 mm WL			
8642.6812	HySafe micro scissor, 5 Fr., 340 mm WL			
Accessories				
823005	BipoTrode electrode, 5 Fr., 360 mm			
8570.040	Reprocessing tray			
8986.251M	Ergonomic handle			
806625301	2.5 mm Light cable kit			
88.01	Sealing cap, 1-2 mm, type S, 12/package			
896.0002	3-Dot stopcock plug, 5/package			



Specifications subject to change without notice.

883-01_05-0618USA