

From the operating room
to the office.



**Diagnostic and Operative
Hysteroscopic Instrumentation**

Bringing care and quality to women's health.

From the OR to the office, our high-precision technologies help you achieve the best clinical outcomes for your patients. Richard Wolf's gynecologic endoscopy portfolio features trusted solutions, such as the 3.9 mm Compact Hysteroscope, HySafe micro-modular hand instruments, and LED 50 HD II Office Camera System, streamlining the transition to office-based procedures. Our advanced instruments are setting new benchmarks in patient safety and ease-of-use.

Optiflow

Continuous-Flow Operating Sheath

Optimized flow dynamics in a single piece design.

The OptiFLOW 5.5 mm Continuous-Flow Operating Sheath features:

- An oblique tip that allows for atraumatic cervical access
- A continuous flow sheath in a single piece design, for improved visualization
- A small footprint, ideal for operating room or office procedures
- Easy assembly and disassembly with 2.7 mm PANOVIEW PLUS telescope



Developed with increased flow dynamics compared to the traditional dual sheath system, the OptiFLOW Continuous-Flow Operating Sheath features a single piece design with a quick-connect lock for simple assembly of the hysteroscope. The 45° angle inflow and outflow channels are ergonomically designed for placement of tubing and access to the entire working length of the sheath.


The built-in, offset 5 Fr. working channel accommodates a broad range of instruments and the atraumatic, oblique tip provides easier access, even in stenotic patients.

Compact Hysteroscopes

3.9 mm and 5 mm


True, continuous flow for brilliant visualization.

Designed with a rod lens optics system, our Compact Hysteroscopes deliver brilliant, detailed visualization. These all-in-one, single piece hysteroscopes feature a built-in, 5 Fr. working channel for expanded treatment options as well as separate inflow and outflow channels for optimal visibility.

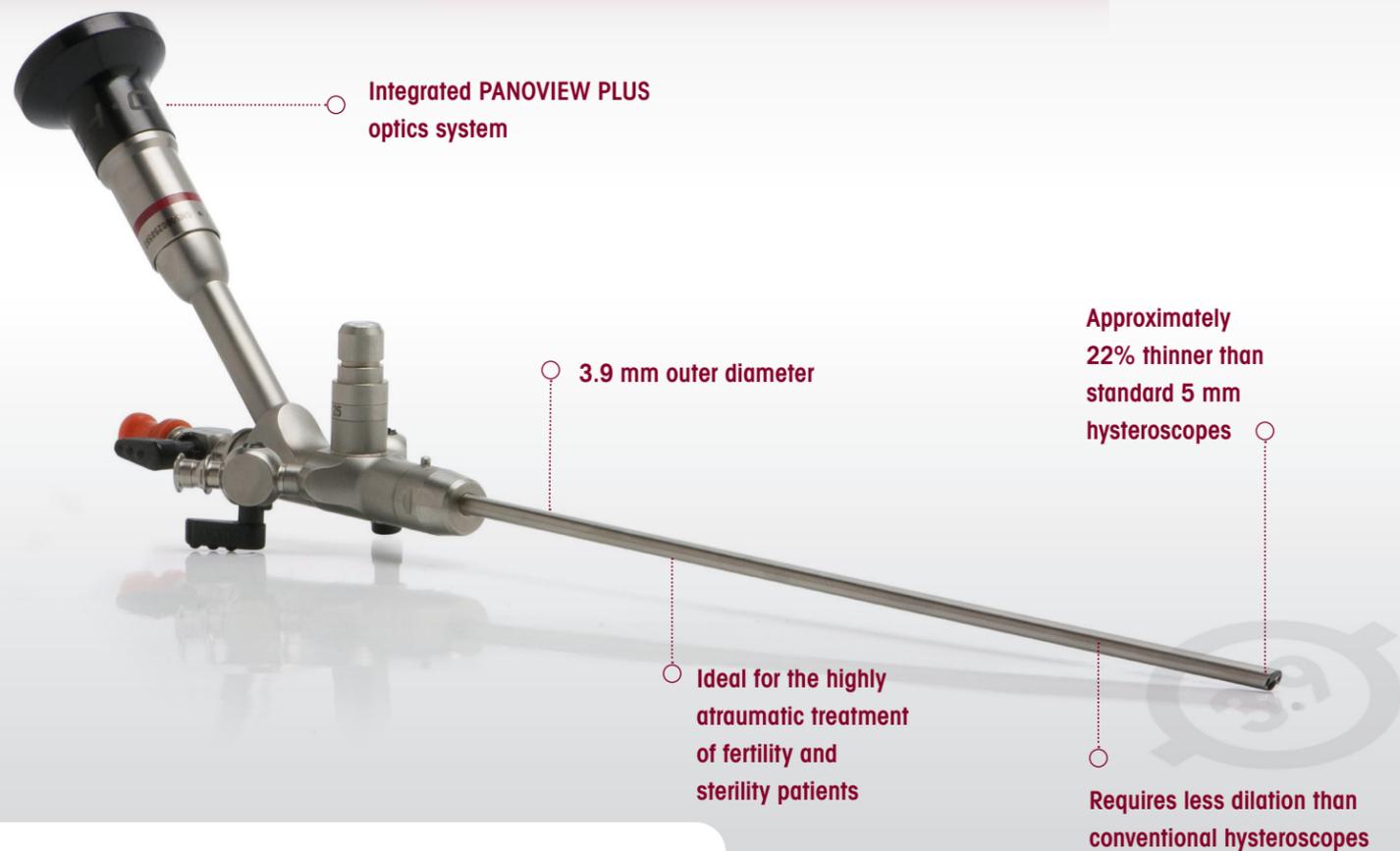

Ideal for use in the office

- With its slim diameter, the 3.9 mm Compact Hysteroscope promotes increased patient comfort and is therefore ideal for use in the office space

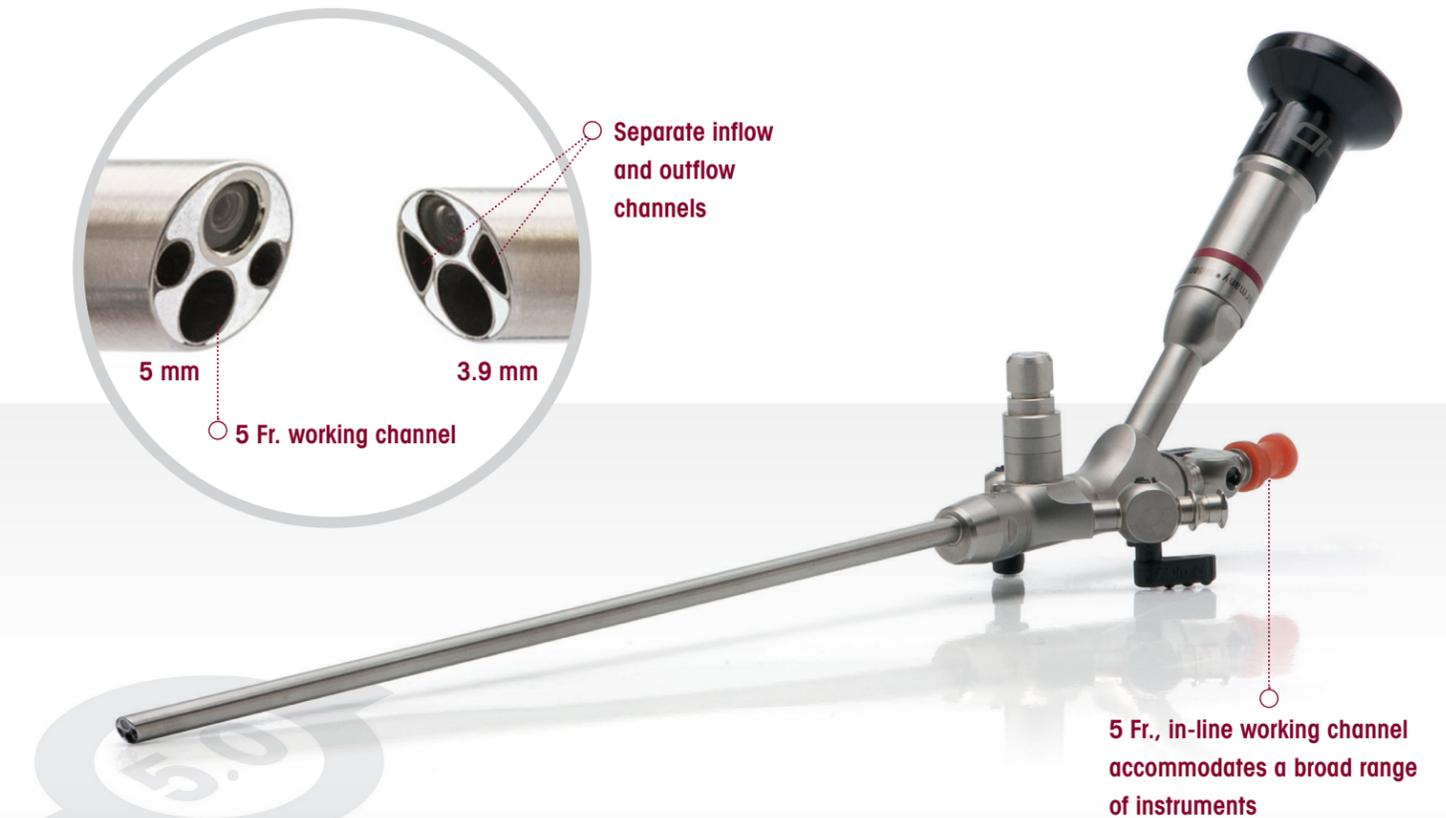

Optional handle for increased ergonomics

- Ergonomic and easy to assemble
- Fits most Richard Wolf hysteroscopes (including the Compact Hysteroscopes and OptiFLOW Continuous-Flow Operating Sheath)



HySafe Micro-Modular Hand Instrument System


360° rotatable, flushable and easy to assemble.



Richard Wolf offers a variety of 5 Fr. hand instrumentation to provide physicians with options for therapeutic intervention during hysteroscopic procedures. The HySafe Micro-Modular Hand Instruments are a 2-piece, modular system that facilitates more thorough reprocessing of individual components and promotes a lower cost of ownership. The 5 Fr. Single Piece Semi-Rigid Hand Instruments require no assembly and offer additional jaw patterns, providing a seamless solution for an array of procedural applications.

5 Fr. Single Piece Semi-Rigid Instruments

Single piece instruments for increased procedural versatility.

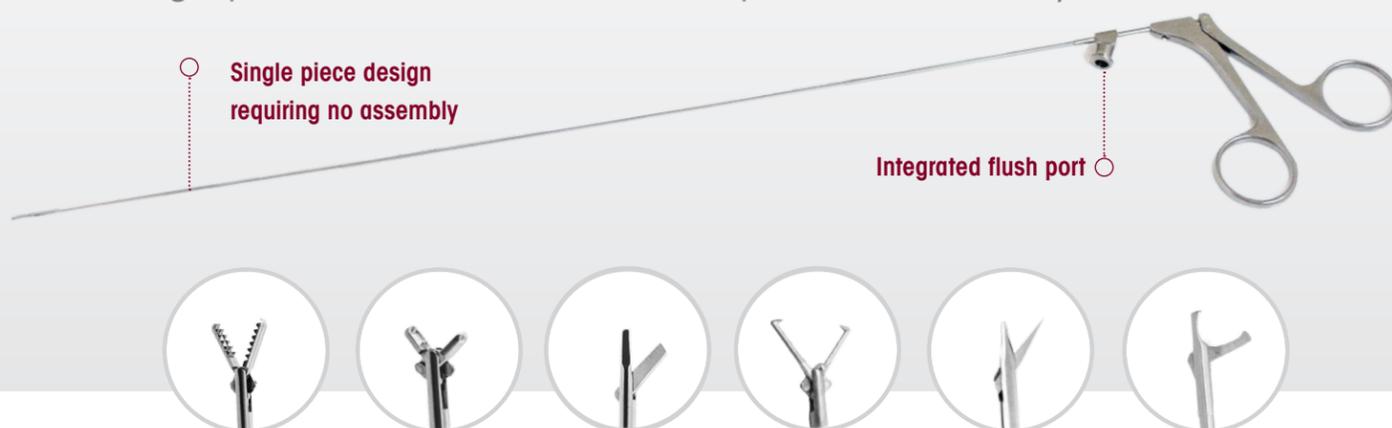

LED 50 HD II Office Camera System

All-in-one.




spirit of excellence

LED 50 HD II All-In-One Unit

- All-in-one unit combining camera controller, light source and integrated recording system allowing for a more economical and space-saving setup optimized for the office
- Plug-and-play for simple setup and intuitive operation
- High quality HD image (1920 x 1080 pixels)
- Bright and efficient LED lamp for homogenous, harmonized illumination
- Integrated recording system for storage of still images and videos on a USB stick

LED 50 HD II Camera Head

- 1/2" CMOS sensor with HD resolution (1920 x 1080 pixels)
- Three programmable camera head buttons for quick and easy access to three functions

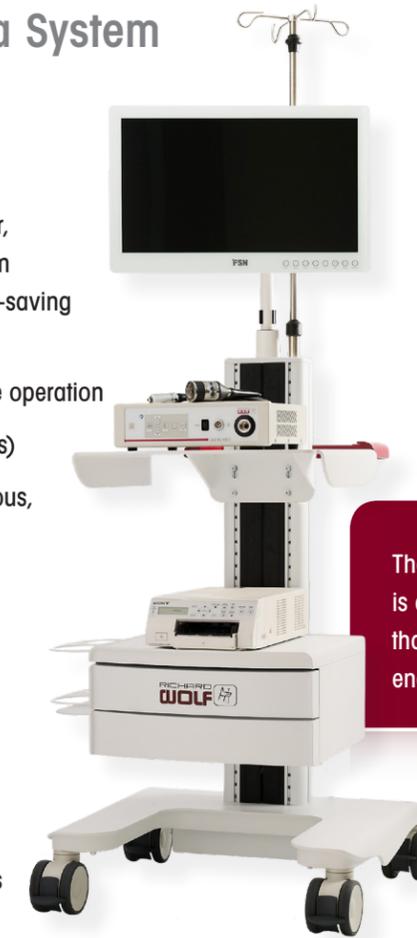

The LED 50 HD II Office Camera System is a space-saving, economical system that provides brilliant, high quality endoscopic imaging.

Office Hysteroscopy Packs

Disposable packs for in-office procedures.

Standard Office Hysteroscopy Pack

(Sterile and disposable – 10 sets per box)

- | | |
|------------------------|---|
| ■ Under buttocks drape | ■ Prep kit:
2 large sponge sticks
2 PVP paint pouches |
| ■ Half drape | |
| ■ Camera drape | ■ Alcohol swab |
| ■ Table cover | ■ Cysto tubing set |
| ■ Blue towels | ■ Extension IV tubing set |
| ■ Leggings | |



Compact Office Hysteroscopy Pack

(Sterile and disposable – 10 sets per box)

- | | |
|------------------------|---------------------------|
| ■ Under buttocks drape | ■ Cysto tubing set |
| ■ Half drape | ■ Extension IV tubing set |
| ■ Camera drape | |
| ■ Table cover | |

Ordering Information

Hysteroscopes & Sheaths

Description	Product No.
OptiFLOW 5.5 mm Hysteroscope Sheath Bundle	89743516
PANOVIEW PLUS Telescope, 12°, Ø 2.7 mm, SL 310 mm	8974.412
PANOVIEW PLUS Telescope, 30°, Ø 2.7 mm, SL 310 mm	8974.402
3.9 mm Compact Hysteroscope Bundle	89864026
5 mm Compact Hysteroscope Bundle	89864016

HySafe Micro-Modular Hand Instruments

Description	Product No.
HySafe 5 Fr. Grasping Forceps Bundle (Alligator)	8642.6502
HySafe 5 Fr. Grasping Forceps Bundle (Spoon Biopsy)	8642.6002
HySafe 5 Fr. Scissors Bundle	8642.6812
HySafe 5 Fr. Biopsy Forceps Bundle (Punch Biopsy)	8642.6312

Single Piece Semi-Rigid Hand Instruments

Description	Product No.
5 Fr. Grasper, Double Action, Single Piece (Alligator)	86422436
5 Fr. Cup Biopsy, Double Action, Single Piece	86422430
5 Fr. Pointed Scissors, Single Action, Single Piece	86422335
5 Fr. Scissors, Single Action, Single Piece (Blunt)	86422434
5 Fr. Hook Scissors, Single Action, Single Piece	86422416
5 Fr. Tenaculum, Double Action, Single Piece	86422341

Office Video System

Description	Product No.
LED 50 HD II All-In-One Unit	862-2519
LED 50 HD II Camera Head	740-2021
RIWO Zoom Lens (f=13-29 mm), Non-autoclavable	5261504
RIWO Fixed Lens (f= 24 mm)	85261244
21.5" HD Medical Monitor	FS-2102D
Monitor Stand	SU-600MD
Condor Multi-Use Equipment Cart	2208030
Condor Cart Video Kit	2208029
UP-27MD Printer	UP-27MD

Recommended Accessories

Description	Product No.
Fiber Light Cable Bundle Ø 2.5 mm, 3 m long	806625301
Ergonomic Handle for Hysteroscope	8986.251
Hysteroscopy Pack (Standard)	4974.0010
Hysteroscopy Pack (Compact)	4975.0010

To arrange a trial, contact your local Richard Wolf Sales Representative or call Customer Service at (800) 323-9653.

Specifications subject to change without notice.

1893-01.03-1225USA

GUARDIAN
Protecting Performance



SERVICE BY RICHARD WOLF