

Size matters



PRINCESS

Slim 7 mm resectoscope
for minimal dilatation effort

PRINCESS

Slim 7 mm resectoscope for minimal dilatation effort

The **Princess** is the uniquely petite resectoscope from R. Wolf with an outer diameter of just 7 mm. It is used with the 12° PANOVIEW Telescope. Very fine and exact resection is possible with these small dimensions, the dilatation effort and the insertion resistance are significantly reduced. The intrauterine pressure can be kept particularly low due to the reduced space requirement. Less fluid absorption and reduced fluid loss are the result.

The petite dimensions of the **Princess** make it ideal for use in

- Stenotic cervical canal
- Fertility and sterility patients



Double sheath system

permits 360° rotation and ensures more freedom of movement.




Classic resectoscope



7 mm Princess Resectoscope

Size matters


- Reduced dilatation effort
- Significantly less insertion resistance
- Substantially lower traumatization to the cervical canal

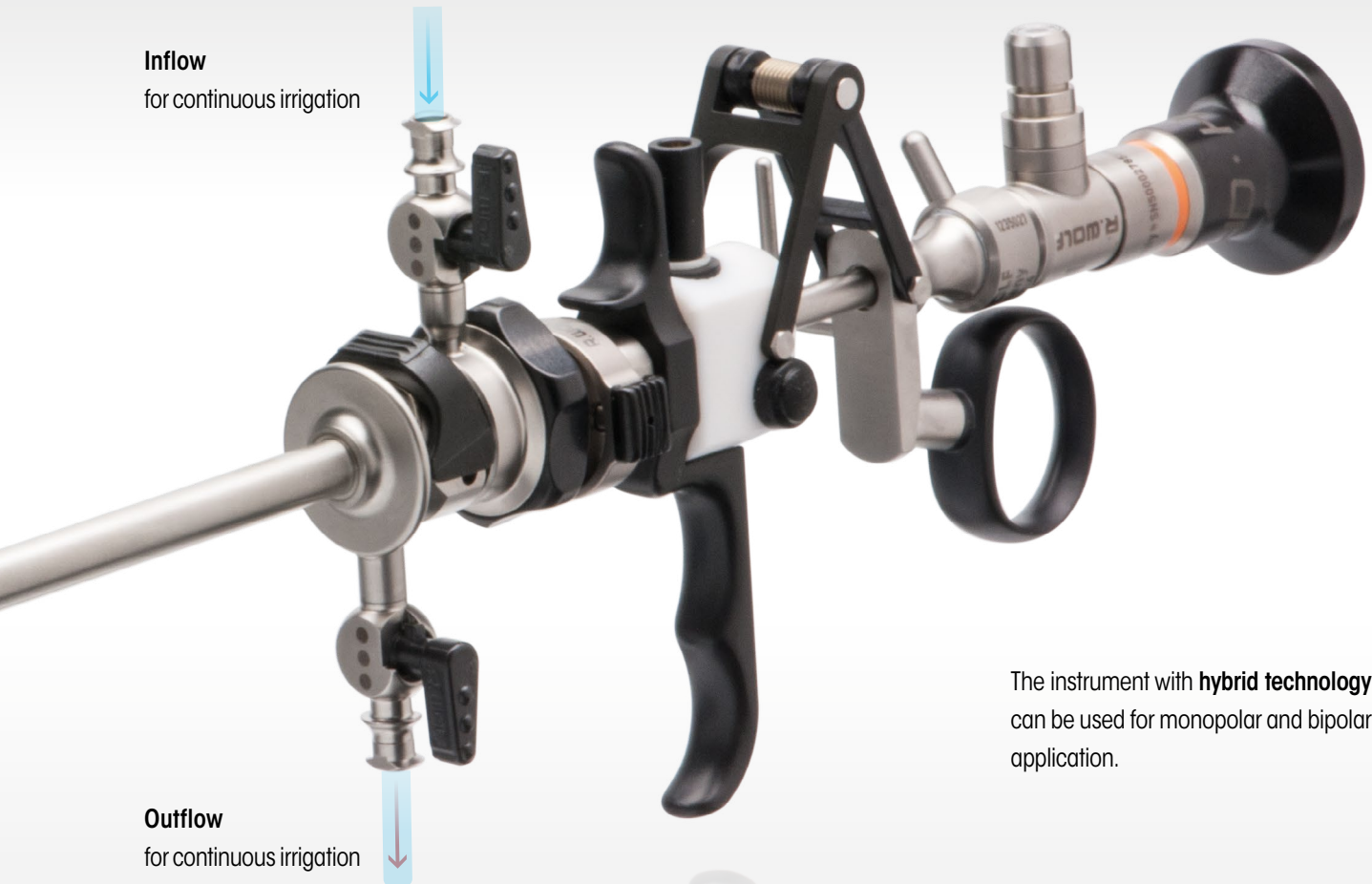

The instrument with **hybrid technology** can be used for monopolar and bipolar application.

Bipolar vaporization electrode "BiVAP"

The new BiVAP bipolar vaporization electrode is ideal for endometrium ablation and treatment of myomas. The 7 mm **Princess** System is also impressive with the significantly reduced space requirement.

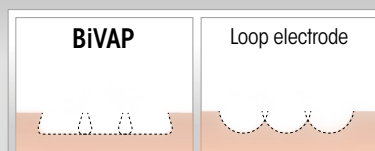

- Optimized electrode shape
- Incremental contact surface for low energy density



- Fast and effective tissue vaporization
- Plasma is concentrated by the unique shape

Size	Version	Outer sheath	Inner sheath	Working element with open handle		Telescope ø 2.7 mm 12°	Obturator	
				passive cutting	active cutting		Viewing obturator	Obturator
7 mm	rigid	8653.021	8653.018	8653.224	8653204	8974.412	865319	8653.18
	rotatable	8653022	8653019					
8 mm	rotatable	8655.334	8655.344					8654.16

Electrodes

		Multi-use	7 mm	8 mm	Types
 Cutting electrode	bipolar, pack of 3, packed singly, sterile		•		4653.1313
	monopolar	•	•		8653.131
 Cutting electrode	bipolar, pack of 3, packed singly, sterile			•	4653.1323
	monopolar	•		•	8653.132
 Roller electrode	bipolar, pack of 3, packed singly, sterile		•	•	4653.023
	monopolar	•	•	•	8653.02
 Hook electrode	monopolar and bipolar	•	•	•	865303
 Hook electrode	monopolar and bipolar, pack of 3, packed singly, sterile		•	•	4653033
 BiVAP electrode	bipolar, pack of 3, packed singly, sterile		•		4653.0213
				•	4653.0243

Optional accessories

Working insert for auxiliary instruments
up to 7 Fr. 8653284

HF bipolar connection cable
for Richard Wolf HF-GENERATOR 2260,
3 m..... 8108233
5 m..... 8108253
for Martin,
3 m.....8108.231
5 m 8108.251
for Erbe,
VIO 300D 20196-119
VIO 3 21196-119

Flexible fiber light cable
Fiber bundle Ø 2.5 mm
(color marking green),
length 2.3 m 806625231

Sterilization protection sleeve
for electrodes 33030

Container (W x H x D)
400 x 60 x 200 mm
incl. lid and peg-mat 33022